



WITH AGILENT, YOU HAVE THE FLEXIBILITY TO CHOOSE THE AUTOMATIC SAMPLER PERFECTLY SUITED FOR THE NEEDS OF YOUR LAB

Whether you're detecting emerging contaminants in drinking water, testing drugs for purity, or testing for illegal drugs, your results have to be correct, precise, and irrefutable. As your lab's requirements and needs continue to change, achieving that level of precision can be challenging.


Agilent can help you meet that challenge head on. With our portfolio of Sample Introduction products, we have a solution that's suited to the specific types of samples and analyses you perform. Regardless of your sample introduction needs, Agilent has a sampler to support your lab's productivity.

Our most flexible liquid sampling system

Agilent 7693A Series Automatic Liquid Sampler (ALS)

Using today's latest technology, the Agilent 7693A helps to ensure your best chromatography with advanced injection features. Whether you have a few or hundreds of samples to analyze, you'll experience sample handling and injection capabilities that are best in class.


- · Multiple-wash solvent capability for pre- and post-injection needle rinsing reduces sample carrvover
- Pre-injection sample pumps and washes further reduce the possibility of carryover
- Single stroke injection volume, as small as 0.01 µL and as large as 250 µL, lets you precisely match injections to the exact needs of your analyses

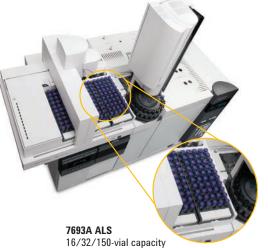







sampling

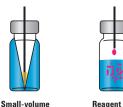

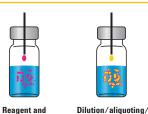

reconstitution



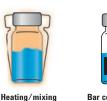

Extended sampling flexibility and automation options. Variable needle depths allow you to sample anywhere within the vial. while optional sample handling capabilities further enhance your productivity.

First-class precision, reliability, and ease of use

Agilent 7697A GC Headspace Sampler

With best-in-class technology and powerful software, the Agilent 7697A Headspace Sampler is packed with the latest productivity-boosting features.

standard addition

- Unique sampling design allows you to use hydrogen as a carrier gas, delivering optimal chromatography and helping future-proof your lab
- · Comprehensive software goes beyond sample handling to guide you through tasks such as method development and resource conservation
- · Method optimization tools facilitate headspace method development
- Electronic pneumatic control, vial leak checking, and barometric pressure compensation ensure consistent results



7697A Headspace Sampler 12/111-vial capacity

Robust productivity tool for intermediate capacity needs **Agilent 7650A ALS**

For labs with intermediate sampling needs, the Agilent 7650A is the ideal solution for achieving maximum productivity.

- State-of-the-art sample handling provides the highest level of reliability for GC sampling
- Delivers performance consistent with other Agilent ALS products
- Programmable injection means plunger speed can be optimized for large volume or challenging applications

Designed specifically for your 6850 GC

Agilent 6850 ALS

Combine the Agilent 6850 Automatic Liquid Sampler (ALS) with the Agilent 6850 GC to create the ideal gas chromatography system for routine QA/QC analyses.

- Programmable, allowing unattended operation for increased productivity
- · Preserves the GC system's small footprint, while adding significant sample capacity
- Optional turret to accommodate 22 4-mL vials

7650A ALS 50-vial capacity 6850 ALS 27-vial capacity

Agilent PAL Sampler

- Combines liquid, headspace, and SPME injection into one instrument, so you can switch from one application to another on the same GC workstation
- Accommodates 4 vial sizes: 1 mL, 2 mL, 10 mL, and 20 mL

Thermal Desorbers

 Labor-saving sample introduction technique for GC and GC/MS that combines selective concentration enhancement with direct extraction into the carrier gas

Tekmar Purge & Trap

- Ability to vary flow rates in different modes, which increases sample throughput
- Ultra-low carryover; superior water management



Agilent PAL Sampler 80 294-vial (2 mL) and 96-vial (10/20 mL) capacity

Agilent PAL Sampler 120 686-vial (2 mL) and 224-vial (10/20 mL) capacity



7667A Mini TDTube-only one-stage technology



Markes TD100 Cryogen-free re-focusing with 100 sample capacity



Tekmar Atomx 80-vial capacity

AT A GLANCE: AGILENT SAMPLE INTRODUCTION DEVICES MEET YOUR NEEDS

This chart provides a quick reference for identifying which Agilent Sample Introduction solution will complement your lab's current applications.

GC and GC/MS	7890 GC	7820A GC	6850 Series II	6890 GC	5975T	6820 GC
7693A ALS	•	•*	**	•	**	
7697A Headspace Sampler	•	•	•	•	•	
7650A ALS		•				
6850 ALS			•			
PAL Samplers	•		•	•		
Markes TD	•	•	•	•	•	•
7667A Mini TD	***	•	***	***		***
Tekmar Purge and Trap	•	•	•	•	•	•
Capillary Trap Sampler	•				•	
Thermal Separation Probe	•				•	

^{*}Injector only, unless used with 5975E GC/MS

^{***} Remote control mode only



Capillary Trap Sampler: This portable, handheld sampler concentrates trace-level airborne compounds in minutes, allowing you to reduce your sampling speed to just minutes — or seconds — per sample.

Thermal Separation Probe: This TSP provides fast analysis of solid, liquid, and slurry samples with limited or zero sample preparation required.

For more information

Whatever your sample introduction needs, we have an instrument that can maximize your lab's productivity.

Visit www.agilent.com/chem/GC

Buy online: www.agilent.com/chem/store

Find an Agilent customer center in your country: www.agilent.com/chem/contactus

U.S. and Canada 1-800-227-9770 agilent_inquiries@agilent.com

Europe

info_agilent@agilent.com

Asia Pacific

inquiry_lsca@agilent.com

This information is subject to change without notice.

© Agilent Technologies, Inc. 2013 Printed in the USA March 12, 2013 5991-1287EN



^{**} Injector only